

medi PT control

Soft brace for patella tracking control with flexion/extension limitation, made of medi Airtex-material

Indications

- patello femoral pain syndrome
- patella chondroplasty
- retropatella osteoarthritis
- post-operative and post-traumatic irritations
- after lateral release surgery to limit the extension and flexion of the joint and to stabilise the patella
- patella and knee instabilities

Mode of action

- collateral stabilisation of the knee joint
- prevention of hyperextension and limitation of range of motion
- correction of lateral patella mal-tracking
- soft brace with compressive and proprioceptive effect
- patella stabilisation buttress that prevents lateral movement throughout full flexion and extension

Features

- polycentric physioglide hinges based on the principle of 4-chain linkage
- made of breathable and moisture wicking medi Airtex material
- extension limitation at: 0°, 10°, 20°, 30°, 45°
flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90° and 120°
immobilisation at: 0°, 10°, 20°, 30°, 45°
- exact correction of the patella-glide due to an individually adjustable patella stabilisation buttress and air-cushion
- separately available Airtex sleeve




stabilisation buttress




correction of patella-glide

Product benefits

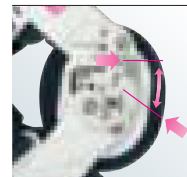


physioglide® hinge

A unique polycentric hinge design that closely mimics the roll and glide motion of the knee.

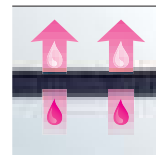

medi airtex

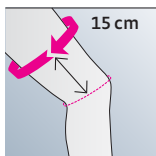

medi PT control

size	circumference cm	article number	
		left	right
XS	31 – 37	G.140.201	G.140.301
S	37 – 45	G.140.202	G.140.302
M	45 – 51	G.140.203	G.140.303
L	51 – 57	G.140.204	G.140.304
XL	57 – 65	G.140.205	G.140.305
XXL	65 – 74	G.140.206	G.140.306

measurement: circumference of thigh 15 cm above the centre of the patella (measured at the medial side of the thigh)