

Carefree Premier

Static Mattress Replacement



The Carefree Premier Mattress is suitable for users who are considered to be up to a 'very high risk' of pressure ulcer development.

The mattress comprises of a contoured "keyhole cut" CMHR foam insert. This castellated foam section provides the best possible pressure reduction together with offering excellent patient comfort and product longevity.


The u-shaped base and side wall construction offer stability and support for patient transfers. The mattress cover has high frequency welded seams and a fully concealed zip which reduces the possibility of fluid ingress thus aiding Infection Control protocols



CMHR foam insert



Reinforced CMHR walls

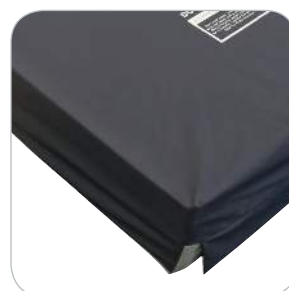


PU cover with full-length flap



Contoured CMHR foam

Features	Benefits
CMHR foam insert	Offers excellent pressure reduction, comfort and support.
Reinforced base and walls	Offers support and assists in patient transfers.
Welded cover	Eliminates fluid ingress and supports infection control protocols.
Toughened PU base	Provides a robust protective surface against hydrolysis / strike through.
2-way stretch, waterproof, vapour permeable PU cover	Improves patient comfort and maintains skin integrity by reducing shear and friction forces whilst aiding infection control.
Concealed zip	Full-length concealed zip prevents fluid ingress and contamination of foam core - to comply with infection control standards.
Non-turn mattress	Does not require turning or rotating, significantly minimising manual handling issues.
Flexible size options	Variable sizes are available upon request.

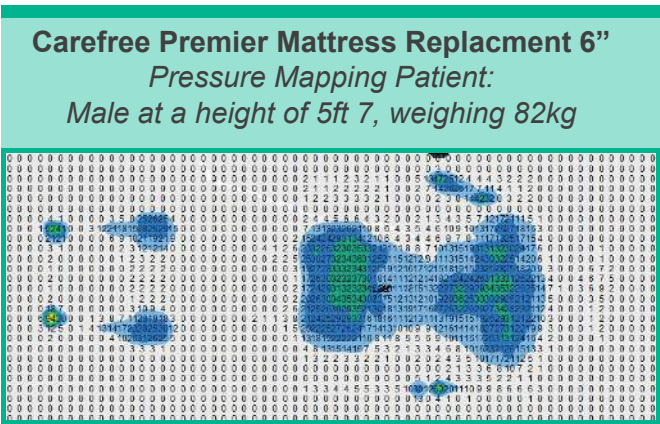

Technical Specification

Product code:	UPRFP3678	Product code:	UPRFP5478-4
Risk category:	Up to very high risk	Risk category:	Up to very high risk
Max user weight:	260 kg / 41 stone	Max user weight:	260 kg / 41 stone
Warranty:	2 Years	Warranty:	2 Years
Mattress dimensions:	1980mm x 880mm x 150mm	Mattress dimensions:	1980mm x 910mm x 100mm
Product weight:	13kg	Product weight:	19kg
Foam density:	Blue: 38-40 kg/m3 White: 38-40 kg/m3		
Nominal hardness:	Blue: 200 Newtons White: 120 Newtons		
Fatigue class:	Very Severe		
Fire retardancy:	BS7177 (Crib 5 Composite)		
Cleaning:	Can withstand up to 10,000 ppm chlorine (to be used in conjunction with local infection control guidelines) (can be laundered up to 84°C)		

Ultimate Healthcare,
 Calmore Industrial Estate, Nutsey Lane, Totton, Hampshire, SO40 3WW
 tel: 0333 321 8996, fax: 0845 051 0035, email: sales@ultimatehealthcare.co.uk
www.ultimatehealthcare.co.uk

