



# COMVI™ Tracheobronchial stent



## SALES OFFICE



Bldg.C 4F, Daebang Triplaon, Jungsan-dong, Ilsandong-gu, Goyang-si, Gyeonggi-do, Korea 410-315  
TEL : +82.31.904.6153 FAX : +82.31.904.6157 E-MAIL : sales@stent.net www.stent.net



# Taewoong Niti-S™ Comvi™ Tracheobronchial stent

## Features

- **Structure** : Triple layered construction(wire-membrane-wire)
  - provide an excellent radial force to maintain luminal patency
  - outer wire mesh prevents the risk of migration
  - biocompatible polytetrafluoroethylene(PTFE) membrane prevents tumor invasion
  - retrieval string at the proximal end : reposition or removal of stent
- **Radiopaque marker** : 4(four) at both ends & 2(two) at the middle

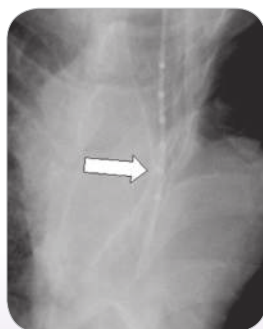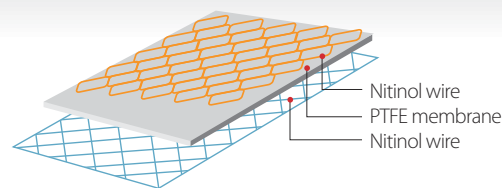

Before stenting : narrowing of the trachea



After stenting

Treatment of Tracheobronchial Obstruction with a Covered stent - Ji Hoon Shin, MD (J Vasc Interv Radiol 2006 ; 17:657-663)

## Ordering Information

### OTW(Over The Wire)

#### ▶ Tracheal stent (Fully covered)

Code	Stent		Delivery				
	Diameter (mm)	Length (cm)	Profile (Fr)	Usable Length (cm)			
TR1603FP	16	3	16				
TR1604FP		4					
TR1605FP		5					
TR1606FP		6					
TR1607FP		7					
TR1608FP		8					
TR1609FP		9					
TR1610FP		10					
TR1803FP		18			3	18	
TR1804FP					4		
TR1805FP	5						
TR1806FP	6						
TR1807FP	7						
TR1808FP	8						
TR1809FP	9						
TR1810FP	10						
TR2003FP	20		3	18	60		
TR2004FP			4				
TR2005FP		5					
TR2006FP		6					
TR2007FP		7					
TR2008FP		8					
TR2009FP		9					
TR2010FP		10					
TR2203FP		22	3			18	
TR2204FP			4				
TR2205FP	5						
TR2206FP	6						
TR2207FP	7						
TR2208FP	8						
TR2209FP	9						
TR2210FP	10						
TR2403FP	24		3	18			
TR2404FP			4				
TR2405FP		5					
TR2406FP		6					
TR2407FP		7					
TR2408FP		8					
TR2409FP		9					
TR2410FP		10					

#### Coding

TR----FP : F-fully covered P-distal side flare

### OTW(Over The Wire)

#### ▶ Bronchial stent (Fully covered)

Code	Stent		Delivery				
	Diameter (mm)	Length (cm)	Profile (Fr)	Usable Length (cm)			
BR1003F	10	3	12	60			
BR1004F		4					
BR1005F		5					
BR1006F		6					
BR1203F		12			3	12	60
BR1204F					4		
BR1205F	5						
BR1206F	6						
BR1403F	14		3	12	60		
BR1404F			4				
BR1405F		5					
BR1406F		6					

#### Coding

BR----F : F-fully covered

#### Relative Studies

- "Benign Tracheobronchial Strictures: Long-Term Results and Factors Affecting Airway Patency After Temporary Stent Placement"
 

- Jin Hyung Kim, et al : AJIR, 2007
- "Placement of Covered Retrievable Expandable Metallic Stents for Pediatric Tracheobronchial Obstruction"
 

- Ho Young Song, et al : JVIR, 2006