


# ELRA™

Endo Luminal Radiofrequency Ablation electrode

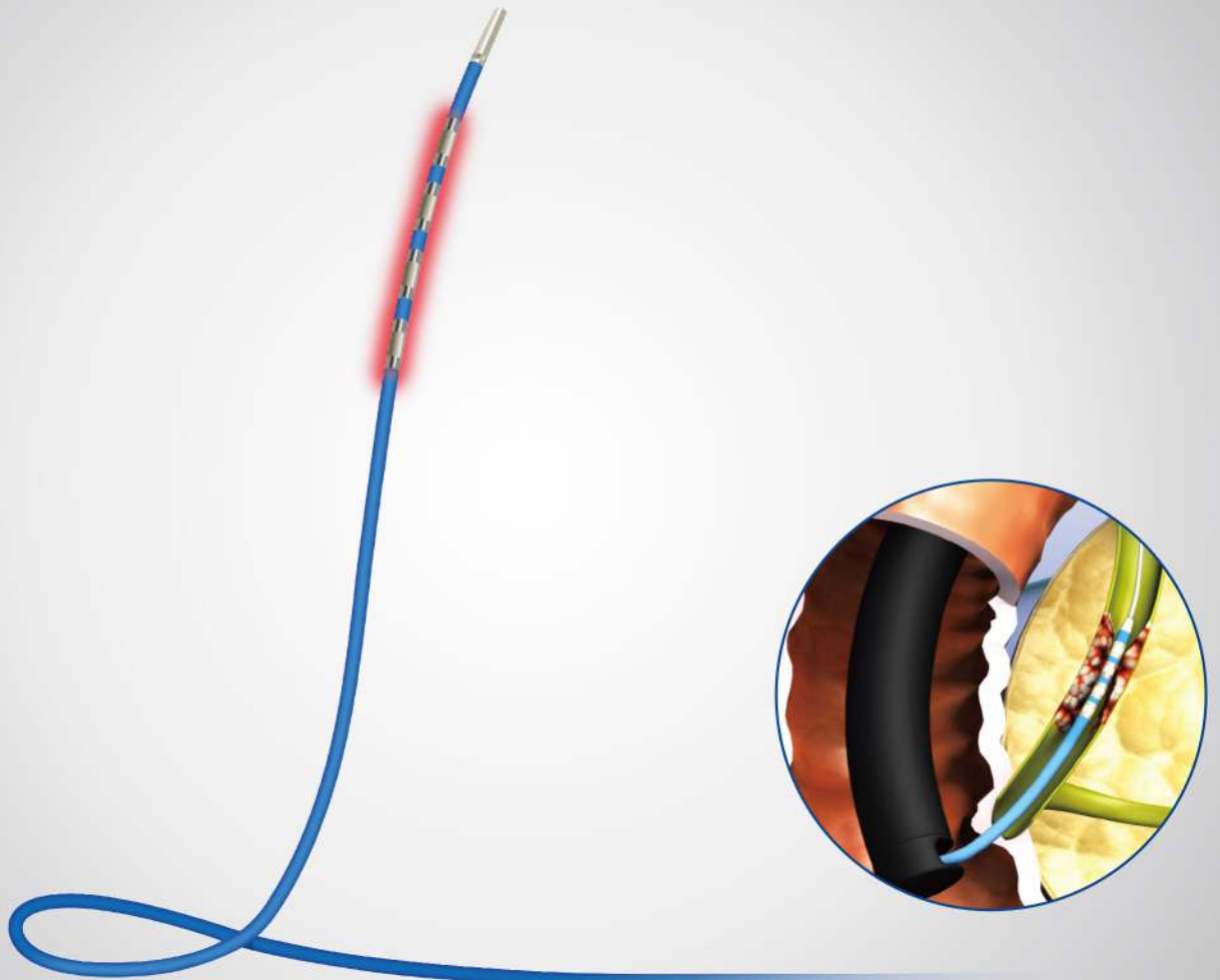

**Global representative** 

Bldg.C 4F, 158 Haneulmaeul-ro, Ilsandong-gu, Goyang-si, Gyeonggi-do, Korea 10355  
E: [contact@stent.net](mailto:contact@stent.net) T: +82 31 904 6153 F: +82 31 904 6157

**Manufacturer** 

Bldg.B 4F, 158 Haneulmaeul-ro, Ilsandong-gu, Goyang-si, Gyeonggi-do, Korea 10355  
E: [info@STARmed4u.com](mailto:info@STARmed4u.com) T: +70 4673 8614 F: +82 0506 816 4546

[www.taewoongmedical.com](http://www.taewoongmedical.com)


**TaeWoong**  
MEDICAL 



for treating malignant biliary obstruction caused by cholangiocarcinoma or pancreatic cancer

## Features


- Optimised temperature based mode makes safe procedure and prevents tissue charring for decreasing any chance of perforation

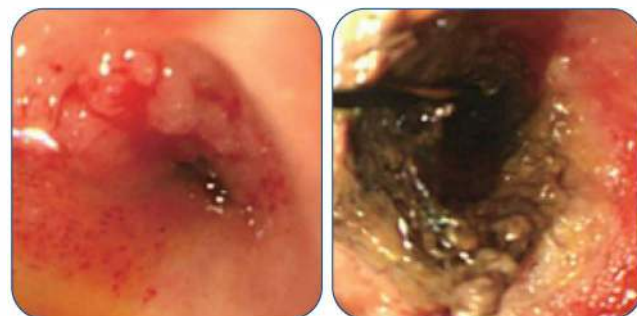


Before Ablation vs After Ablation

### Ablation Image (Bovine Liver)

Setting Condition: 75°C 10watt 2min

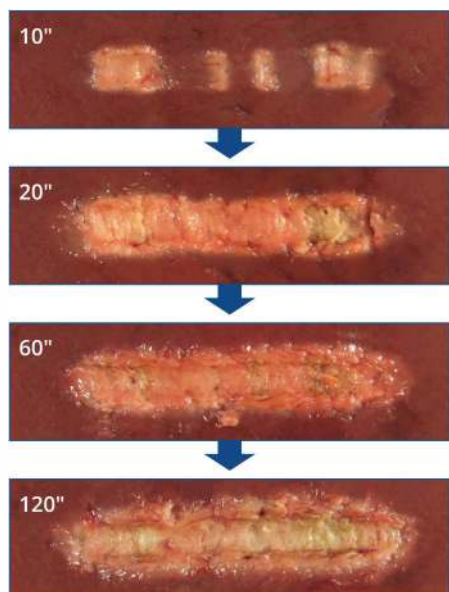

- Multiple bipolar electrode makes linear ablation zone
- Electrode can be selected according to the length of stricture size : 18mm / 33mm

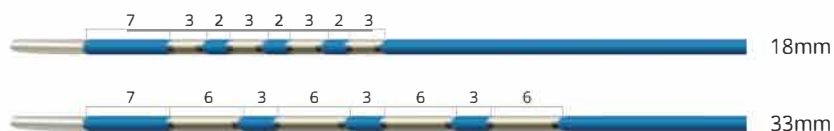

- Electrode is bipolar system and single use ONLY
- Endoscopic and percutaneous approach are available
- It is compatible with VIVA combo™ generator ONLY

## Ordering Information

Endoscopic Approach				Percutaneous Approach			
Code	Total Length (cm)	Exposure Length (mm)	Diameter	Code	Total Length (cm)	Exposure Length (mm)	Diameter
7-4B18L	175	18	7Fr (2.31mm)	7-4B18S	40	18	7Fr (2.31mm)
7-4B33L		33		7-4B33S		33	