


# ELRA™

Endo Luminal Radiofrequency Ablation electrode

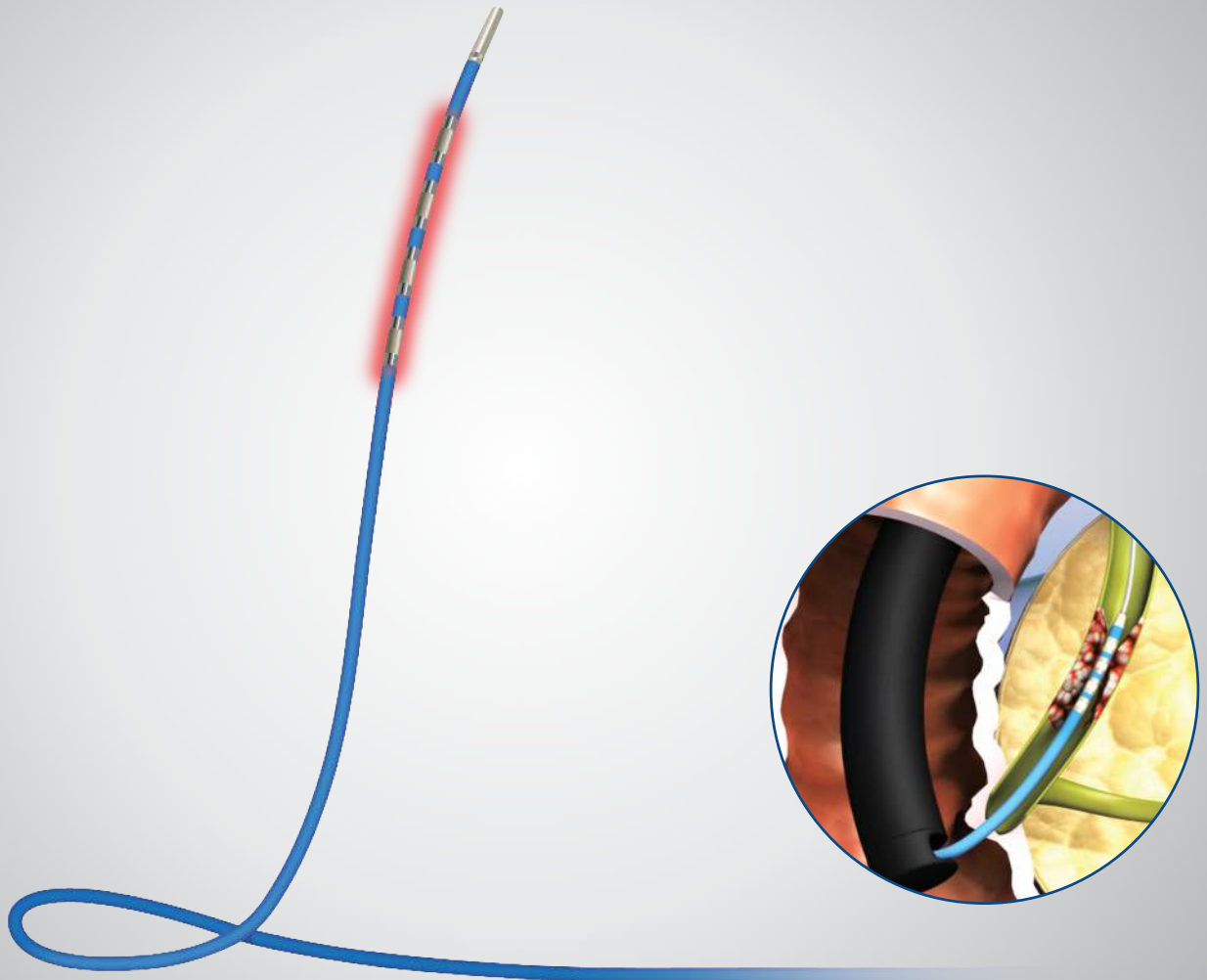

UK Distributor  
Aquilant House, Unit B1-B2 Bond Close,  
Kingsland Business Park, Basingstoke,  
Hampshire, RG24 8PZ, United Kingdom

T: +44 (0) 1256 365 454  
E: [contactus@aquilantservices.com](mailto:contactus@aquilantservices.com)  
W: [www.aquilantservices.com](http://www.aquilantservices.com)

Global representative 

Bldg.C 4F, 158 Haneulmaeul-ro, ilsandong-gu, Goyang-si, Gyeonggi-do, Korea 10355  
E: [contact@stent.net](mailto:contact@stent.net) T: +82 31 904 6153 F: +82 31 904 6157

Manufacturer 

Bldg.B 4F, 158 Haneulmaeul-ro, ilsandong-gu, Goyang-si, Gyeonggi-do, Korea 10355  
E: [info@STARmed4u.com](mailto:info@STARmed4u.com) T: +70 4673 8614 F: +82 0506 816 4546


[www.taewoongmedical.com](http://www.taewoongmedical.com)



for treating malignant biliary obstruction caused by cholangiocarcinoma or pancreatic cancer

## Features


- Optimized temperature based mode makes safe procedure and prevents tissue charring for decreasing any chance of perforation

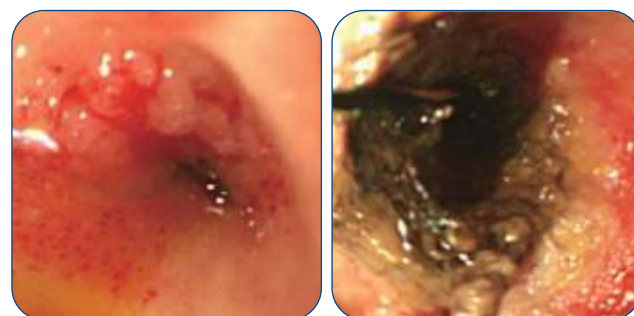


Before Ablation vs After Ablation

### Ablation Image (Bovine Liver)

Setting Condition: 75°C 10watt 2min

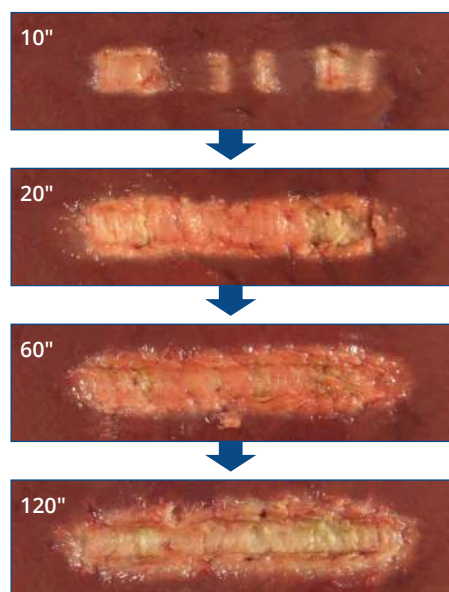

- Multiple bipolar electrode makes linear ablation zone
- Electrode can be selected according to the length of stricture size : 18mm / 33mm

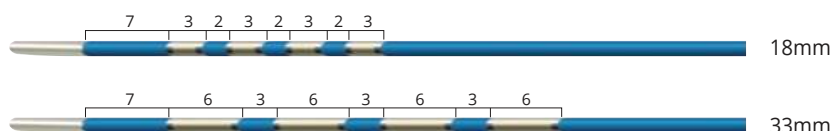

- Electrode is a bipolar system and single use ONLY
- Endoscopic and percutaneous approach are available
- It is compatible with VIVA combo™ generator ONLY

## Ordering Information

Endoscopic Approach			
Code	Total Length (cm)	Exposure Length (mm)	Diameter
7-4B18L	175	18	7Fr (2.31mm)
7-4B33L		33	

Percutaneous Approach			
Code	Total Length (cm)	Exposure Length (mm)	Diameter
7-4B18S	40	18	7Fr (2.31mm)
7-4B33S		33	