

JOBST® Compri2

Two-Layer Bandage System

JOBST® Compri2 is Latex-free and has accessible guidance tools to aid application and promote the correct compression. Two compression levels, ensuring the right level can be selected for both venous and mixed aetiology leg ulcers.

Compression on demand


Two choices to achieve the compression needed

Type of disease	Incidence	ABPI	Compression level	
Venous	~60%	>0.8	40 mmHg	JOBST® Compri2
Mixed aetiology	~20%	0.5 - 0.8	20 - 30 mmHg	JOBST® Compri2 lite
Arterial	~20%	<0.5	No compression	-



1st layer

- Application around the heel becomes much easier because the bandage stretches in both dimensions
- Tearable to allow extra padding to be easily applied where necessary

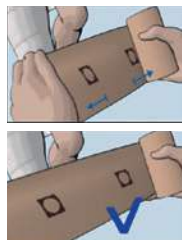

2nd layer

- A simple indicator to assist achieving the correct compression level every time
- The bandage can be adjusted during application



Code	Ankle Size	PIP Code	NHS SC Code
JOBST® Compri2 - regular compression (ABPI >0.8)			
76271-00 (Latex-free)	18 - 25cm	400-6177	EBA145
76271-01 (Latex-free)	25 - 32cm	400-6185	EBA146
JOBST® Compri2 lite - light compression (ABPI 0.5 - 0.8)			
76271-02 (Latex-free)	18 - 25cm	400-6193	EBA147
76271-03 (Latex-free)	25 - 32cm	400-6201	EBA148

