


Novacare High Risk ASX

ASX

Alternating pressure system for the treatment of pressure ulcers up to grade III (according to EPUAP).

The ASX mattress overlay features nylon/PU cells, eight of which are micro-perforated. The cover is bi-elastic, PU-coated and has an all-round zip. In addition, a double layer base permits the insert of a safety mattress of up to 5 cm. This feature protects the existing mattress from soiling.

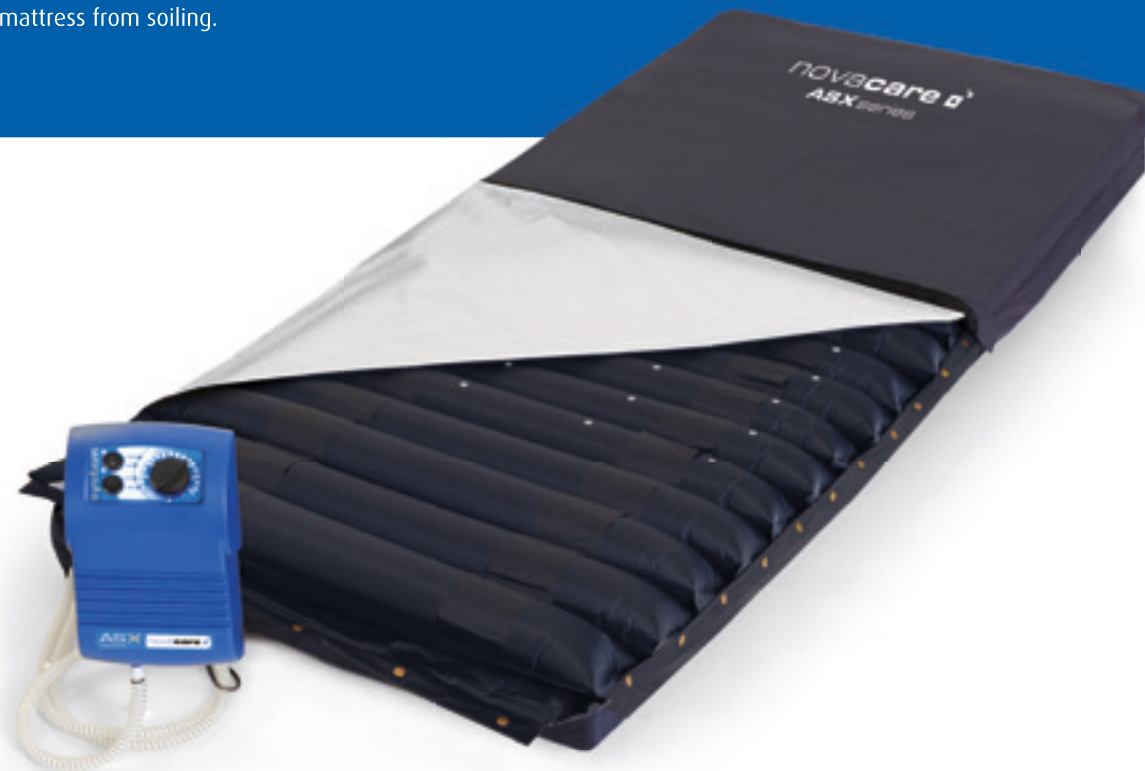

optional insert of mattress base

Mattress Specifications

Item no. 997008

Control unit

Measurements: approx. 30 x 18 x 12 cm
Weight: approx. 2.7 kg

1 membrane compressor, mech. control, alternating pressure cycle (12 min) & static mode, low pressure alarm (visual and acoustic)

Mattress overlay

Measurements inflated (without safety mattress): approx. 200 x 90 x 14 cm
Weight (without safety mattress): approx. 5.4 kg

18 cells, 3 head cells static, non-bending connecting tubes (incl. transportation safety lock), CPR emergency ventilation

Materials

Cells: high-quality nylon PU cell material
Cover: PES polyurethane cover with all-round zip, bi-elastic, highly breathable, slip-retardant

Accessories (optional)

transport bag for mattress and control unit
foam mattress base with polyester cover (item no. 991508+)

Pressure ulcer grade:

I II III IV

Max. patient weight: 140 kg