



Esophageal Dilatation Systems

M-Flex® Blue Silicone Bougies



There's nothing more flexible than an M-Flex®.

Our silicone bougies feature added flexibility that some users prefer. Tungsten-filled and not made with natural rubber latex, there's nothing quite like the M-Flex®.

Benefits

- Ultra-flexible tapered tip (Maloney)
- Completely radiopaque from distal tip to proximal end
- American and European depth markings
- Mercury-free and not made with natural rubber latex
- Autoclavable for quick reprocessing

Features

- 21 sizes: 20 FR to 60 FR (in 2 FR increments)
- Choose a full set, mini set or individual pieces
- Three-year warranty on product and workmanship
- Easy disposal
- No charge trade-in program for proper disposal of expired mercury-filled bougies (with purchase of new product)



Store and transport M-Flex® bougies in our storage case which holds up to two full sets.

M-Flex® Dilatation Systems

Set	Part No.	Contents
Maloney (tapered) 21-piece full set	1206-00	21 bougies, 20-60 FR
Maloney (tapered) 10-piece mini set	1206-10	10 bougies, 36-54 FR
Hurst (blunt) 21-piece full set	1207-00	21 bougies, 20-60 FR
Hurst (blunt) 10-piece mini set	1207-10	10 bougies, 36-54 FR
Bougie Storage Case	1210-01	Heavy-duty storage for up to 46 bougies. 32"L x 7"W x 4"H




Our storage case allows easy identification of bougie size. Removing plexiglass window provides easy access.



End caps display French size.



Extra flexible tip.

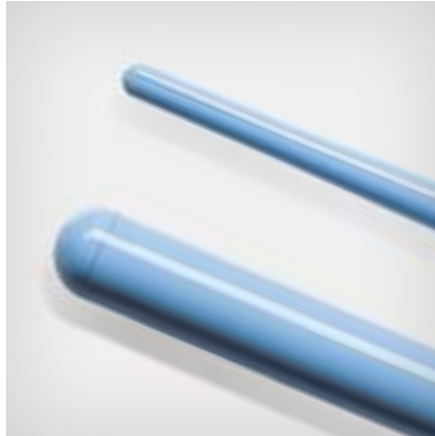

Top: Maloney tip. Bottom: Hurst tip.

Hurst Bougie

Maloney Bougie

American Depth Markings


European Depth Markings



M-Flex® Blue Silicone Bougies

FR (French) Size	(mm)	Maloney (tapered) Tip Part No.	Hurst (blunt) Tip Part No.
20	6.6	1206-20	1207-20
22	7.4	1206-22	1207-22
24	8	1206-24	1207-24
26	8.6	1206-26	1207-26
28	9.3	1206-28	1207-28
30	10	1206-30	1207-30
32	10.7	1206-32	1207-32
34	11.2	1206-34	1207-34
36	12	1206-36	1207-36
38	12.5	1206-38	1207-38
40	13	1206-40	1207-40
42	14	1206-42	1207-42
44	14.5	1206-44	1207-44
46	15.2	1206-46	1207-46
48	16	1206-48	1207-48
50	16.5	1206-50	1207-50
52	17.2	1206-52	1207-52
54	18	1206-54	1207-54
56	18.5	1206-56	1207-56
58	19.2	1206-58	1207-58
60	20	1206-60	1207-60

WeightRight® Mercury-Free Bougies

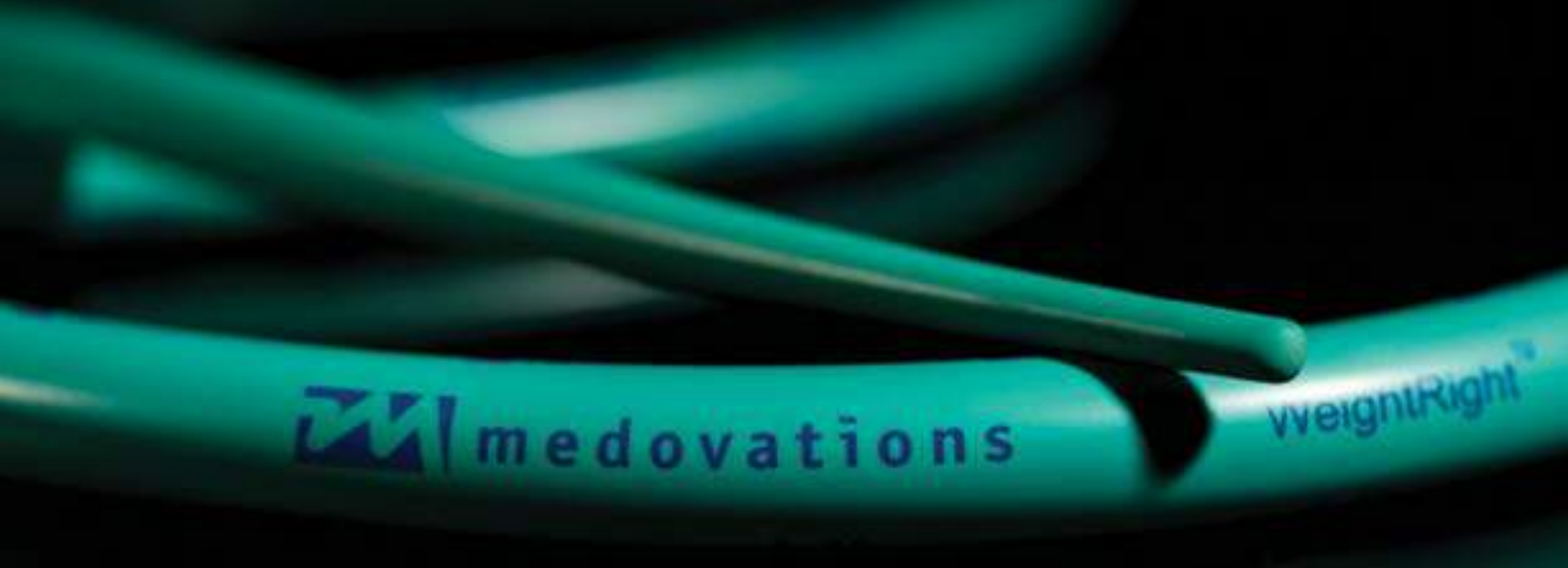

Heavy-weight control and convenience combined.

Our tungsten-filled bougie mimics the weight and feel of mercury-filled bougies without the concerns of mercury usage and disposal. It's the weight you expect, and the convenience you demand.

Benefits

- Completely radiopaque from distal tip to proximal end
- Flexible material does not attract particulate
- American and European depth markings
- Mercury-free and not made with natural rubber latex

Features

- 23 sizes: 16 FR to 60 FR (in 2 FR increments)
- Choose a full set, mini set or individual pieces
- Three-year warranty on product and workmanship
- Easy disposal
- No charge trade-in program for proper disposal of expired mercury-filled bougies (with purchase of new product)



Store and transport WeightRight® bougies in our storage case.

WeightRight® Dilatation Systems

Set	Part No.	Contents
Maloney (tapered) 23-piece full set	1208-00	23 bougies, 16-60 FR
Maloney (tapered) 10-piece mini set	1208-10	10 bougies, 36-54 FR
Hurst (blunt) 23-piece full set	1209-00	23 bougies, 16-60 FR
Hurst (blunt) 10-piece mini set	1209-10	10 bougies, 36-54 FR
Bougie Storage Case	1210-01	Heavy-duty storage for up to 46 bougies. 32"L x 7"W x 4"H



Our storage case allows easy identification of bougie size. Removing plexiglass window provides easy access.



Hang bougie vertically using the hole in the end cap.



Top: Maloney tip. Bottom: Hurst tip.

Hurst Bougie

Maloney Bougie

American Depth Markings

European Depth Markings



WeightRight® Mercury-free Bougies

FR (French) Size	(mm)	Maloney (tapered) Tip Part No.	Hurst (blunt) Tip Part No.
16	5.3	1208-16	1209-16
18	6	1208-18	1209-18
20	6.6	1208-20	1209-20
22	7.4	1208-22	1209-22
24	8	1208-24	1209-24
26	8.6	1208-26	1209-26
28	9.3	1208-28	1209-28
30	10	1208-30	1209-30
32	10.7	1208-32	1209-32
34	11.2	1208-34	1209-34
36	12	1208-36	1209-36
38	12.5	1208-38	1209-38
40	13	1208-40	1209-40
42	14	1208-42	1209-42
44	14.5	1208-44	1209-44
46	15.2	1208-46	1209-46
48	16	1208-48	1209-48
50	16.5	1208-50	1209-50
52	17.2	1208-52	1209-52
54	18	1208-54	1209-54
56	18.5	1208-56	1209-56
58	19.2	1208-58	1209-58
60	20	1208-60	1209-60

SafeGuide® Over the Guidewire Esophageal Dilators



Flexibility and control for exacting procedures.

At times, there is need for greater accuracy. That's why there's SafeGuide®, a remarkably flexible over the guidewire solution that ensures exact placement of the dilator. Position the guidewire with the assistance of an endoscope and direct the SafeGuide® dilator where you need it.

Benefits

- Completely radiopaque from distal tip to proximal end
- Flexible material does not attract particulate
- American and European depth markings
- Not made with natural rubber latex
- Dilators are reusable

Features

- 16 sizes: 15 FR to 60 FR (in 3 FR increments)
- Choose a full set, mini set or individual pieces
- Easy to clean with optional brushes and adapters
- With proper care, reusable SafeGuide® dilators provide years of service




SafeGuide® Dilatation Systems

Set	Part No.	Contents
16-piece full set	1214-00	16 dilators, 15-60 FR, 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters
7-piece mini set	1214-10	7 dilators (15, 21, 27, 33, 39, 42, 45 FR), 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters
Custom mini set	Contact your sales representative	7 dilators of your choice, 1 storage case, 1 marked spring tip guidewire, 1 reusable cleaning brush, 25 disposable cleaning adapters

Store and transport SafeGuide® dilators in a convenient storage case. A labeled slot holds each dilator secure. The case holds a full set of 16 dilators.



The single-piece construction of the SafeGuide® Spring Tip Guidewire resists breakage and flexes easily.

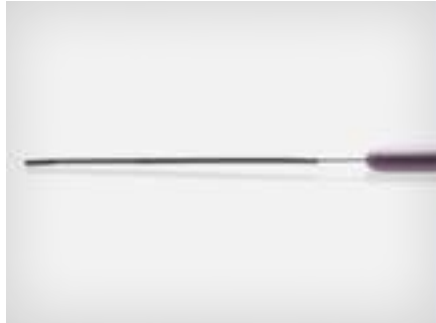

SafeGuide® Over the Guidewire Esophageal Dilators

FR (French) Size	(mm)	Part No.
15	5	1214-15
18	6	1214-18
21	7	1214-21
24	8	1214-24
27	9	1214-27
30	10	1214-30
33	11	1214-33
36	12	1214-36
39	13	1214-39
42	14	1214-42
45	15	1214-45
48	16	1214-48
51	17	1214-51
54	18	1214-54
57	19	1214-57
60	20	1214-60

Accessories



Storage case indicates size.

SafeGuide® Accessories

Set	Part No.	Description	Qty/Case
SafeGuide® Storage Case	1214-01	Holds up to 16 SafeGuide® Dilators. 36"L x 10"W x 4"H	1
Marked Spring Tip Guidewires	1214-02	Guides the dilator. 210cm	2
Disposable Cleaning Brushes	1214-03D	Cleans the inner lumen of the dilator	3
Reusable Cleaning Brushes	1214-03R	Cleans the inner lumen of the dilator	3
Disposable Cleaning Adapters	1214-04	Used to flush the inner lumen of the dilator	25



Clearly marked reusable and single-use cleaning brushes.



Insert cleaning brush into the distal end and push through to the other side for complete inner lumen cleaning.



Disposable cleaning adapter.



Is it time to replace your esophageal bougies or wire-guided dilators?

Replace bougies when you notice:

- End cap separating from bougie
- Cuts, cracks or nicks
- Flaking or peeling of the outer coating
- Faded or missing markings
- Expired or missing expiration dates

Replace wire-guided dilators when you notice:

- Loss of flexibility
- Cuts, cracks or nicks
- Stretching in the inner lumen
- Faded or missing markings
- Curved dilator or curled tapered tip
- Expired



Esophageal Dilators



CleanFreak® Cleaning Products



Endoscope Care & Accessories

At Medovations, we are committed to providing you with high-quality, durable products for your gastrointestinal procedures. See the complete line of GI solutions at Medovations.com

Learn more from your personal product representative at **800-558-6408** or Medovations.com



Medovations, Inc. | 102 East Keefe Avenue | Milwaukee, Wisconsin 53212 USA
800.558.6408 or 414.265.7620 | Fax 414.265.7628
Orders: orders@medovations.com | General Inquiries: medo@medovations.com