


# ULTRAVISION CT Lift IMA Exposure System

---

Leveraging Strength Through Simplicity

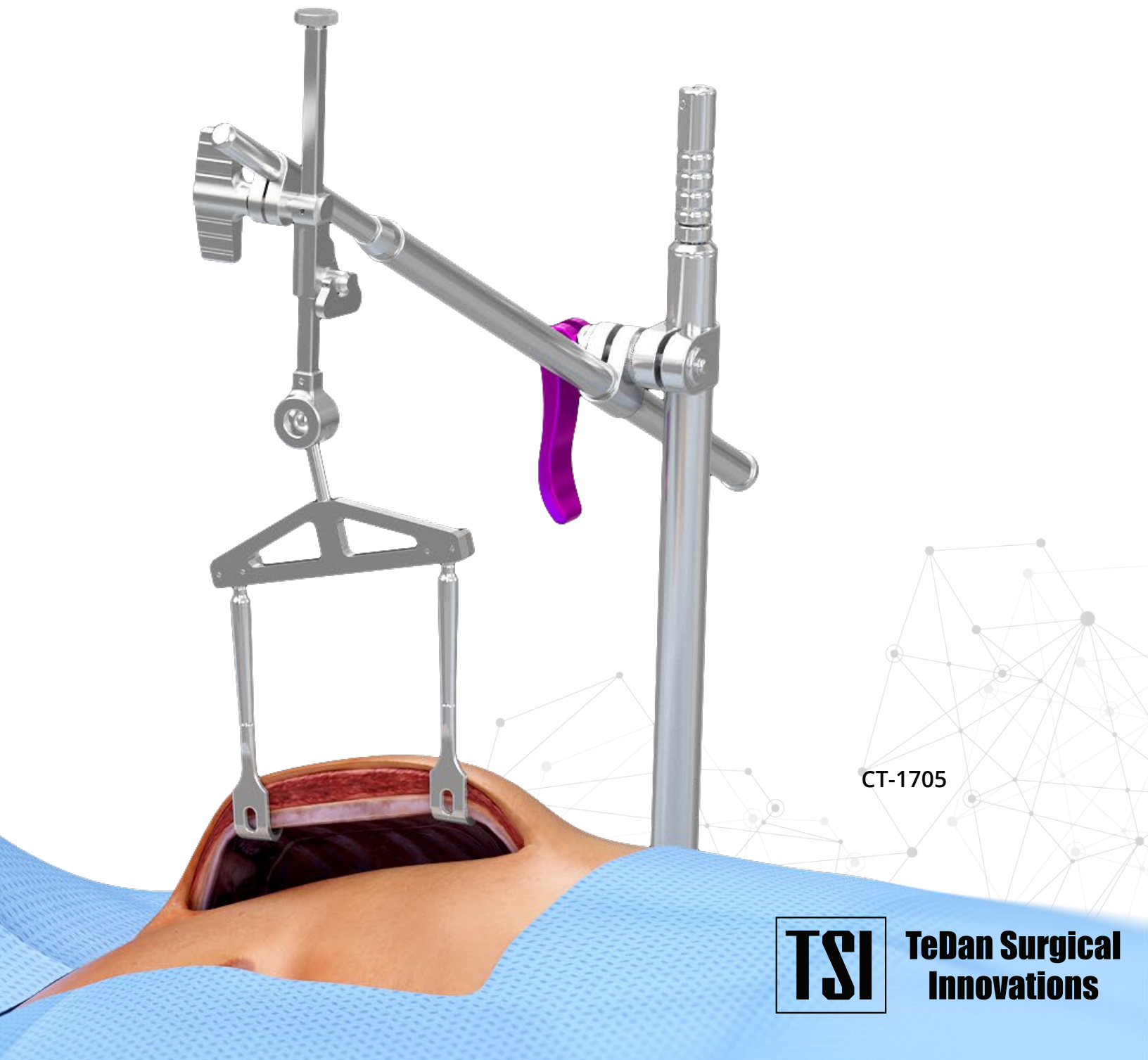

**TeDan Surgical  
Innovations**


CT-1705

ULTRAVISION CT Lift IMA  
Exposure System




# As Simple As 1,2,3


Tedan Surgical Innovations has transformed IMA access for open and minimally invasive CABG procedures with the ULTRAVISION™ CT Lift System. Enhanced with cutting-edge design, it offers optimal exposure and positioning versatility, while also streamlining setup and reprocessing times for faster OR turnaround. The efficient three-step setup minimizes OR staff training requirements, and the system can be set up over the drape, supporting a sterile workflow. The cable-free design eliminates maintenance issues associated with traditional retractors, reducing repair costs and OR downtime.

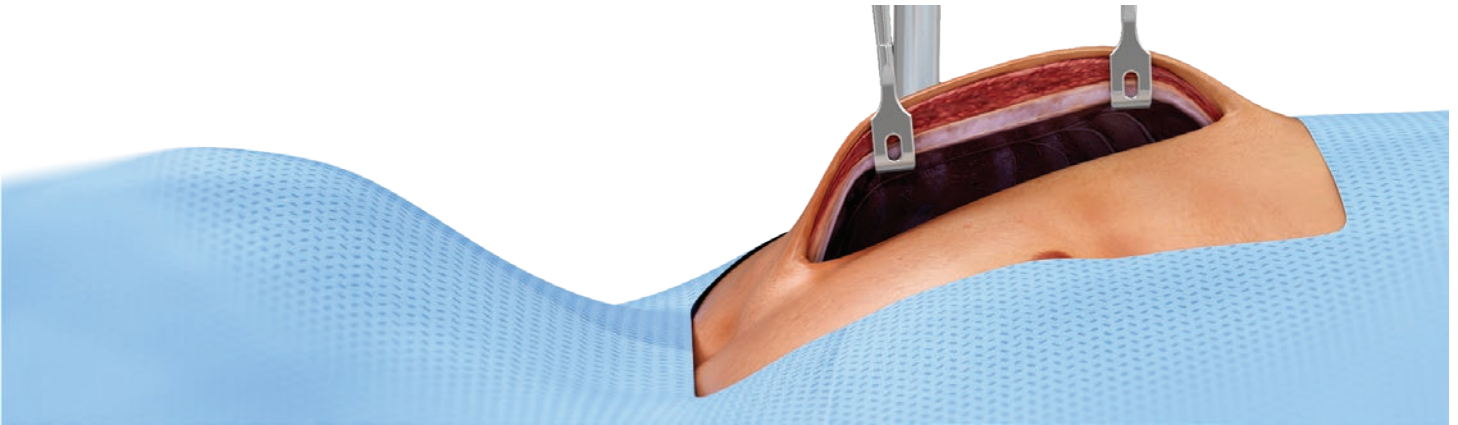

1. Low profile table mounted arm attaches over the drape.

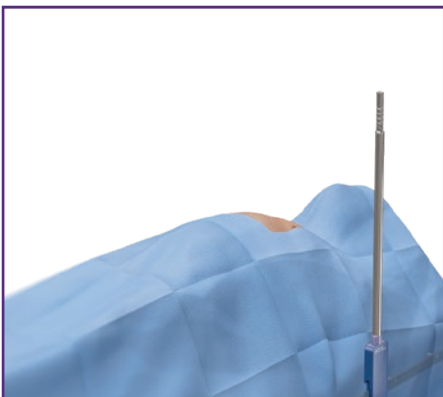

2. Multi directional arm adjustment over patient.

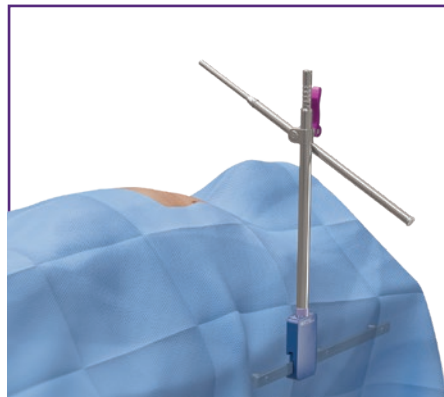

3. Cam-lock mechanism secures set-up in a single step.




# Cable-Free Lift Mechanism

Unlike traditional cable-based systems, the ULTRAVISION CT Lift System employs a rack-based lift mechanism, eliminating issues associated with braided steel cables. This eliminates the risk of frayed cables, reducing the potential for punctured gloves, simplifying cleaning, and minimizing required maintenance.



## Stable and Versatile Design

The system's versatile positioning capability accommodates varying patient anatomy and surgical techniques. Its robust bar design prevents bowing or flexing, ensuring predictable lifting of the chest wall.

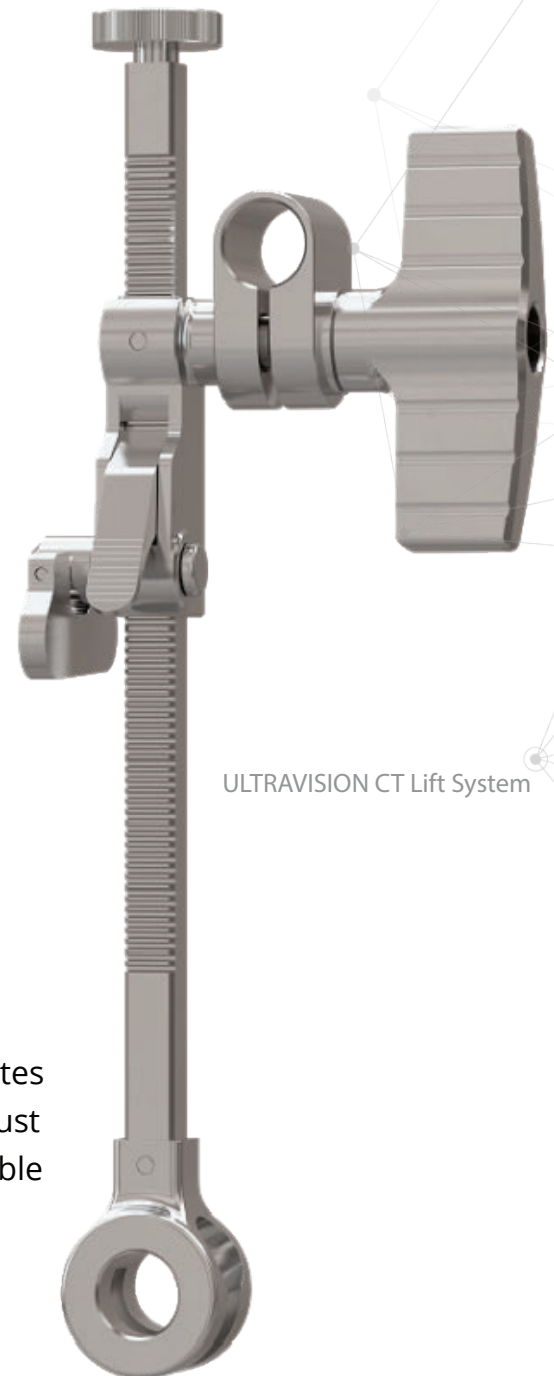



CAT. NO.	DESCRIPTION	QTY
CS-0050	IMA Blade Hook	1
CT-0052	Open IMA Dual Hook	1
CT-0300	IMA Blade, Wide	1
CT-0302	IMA Blade, Narrow	1
CT-0705	Instrument Tray, CT-1705	1
HS-0016	Table Clamp, 450 MM, Purple	1
HS-0019	Extension Arm	1
HS-0100	Lift Assembly, Purple	1


## Accessories

### Optional Items

CAT. NO.	DESCRIPTION	QTY
CT-0105	Universal Stabilizer Adapter	1
CT-0320	Rack Extension Stabilizer Mount	1
HS-0005	Universal Side Rail Extension, 36"	1
LS-0101	T-Handle Crank Control	1

## Illumination

CAT. NO.	DESCRIPTION	QTY
ML-0046	Light Cable, 10', Sterile, 5/BOX	1

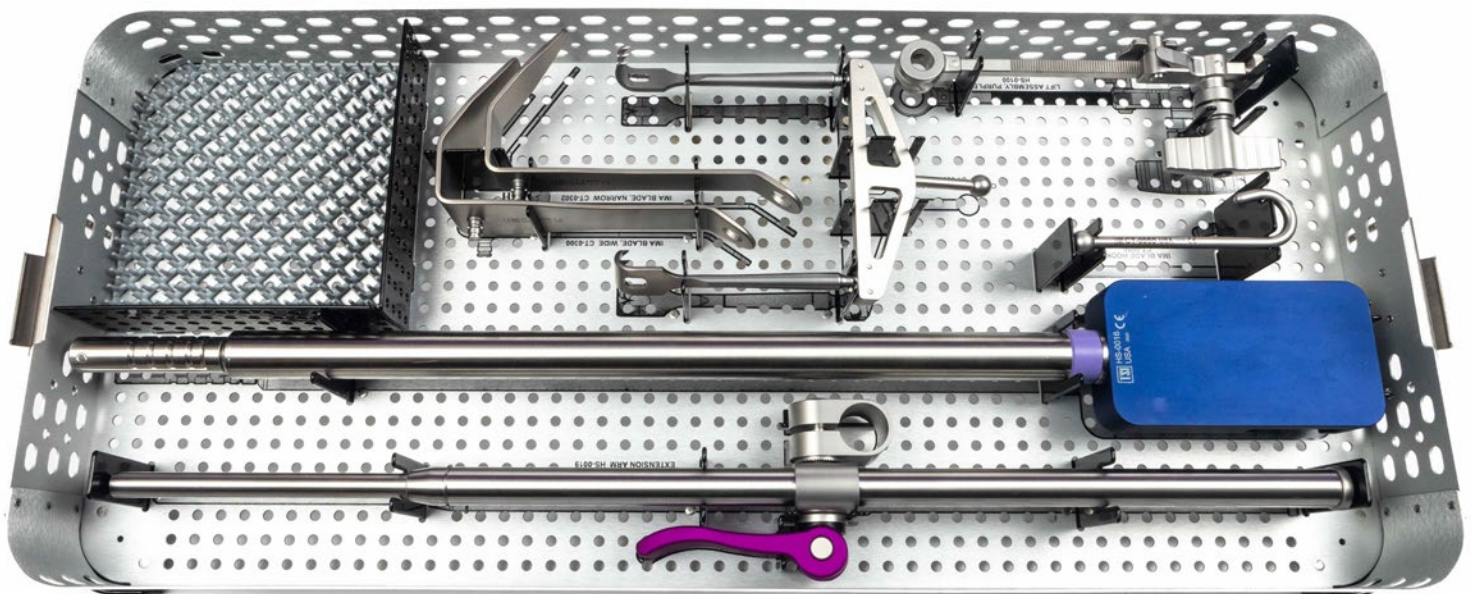



# **TeDan Surgical Innovations**

TeDan Surgical Innovations, Inc.  
12320 Cardinal Meadow Drive, Suite #150  
Sugar Land, TX 77478  
Ph: (+1) 713-726-0886 \* Fax: (+1) 713-726-0846  
Toll Free: 877-726-0886  
[www.tedansurgical.com](http://www.tedansurgical.com)  
[info@tedansurgical.com](mailto:info@tedansurgical.com)

Images for illustration purposes only.

LIT-0089, REV 04 © TSI 2024