



Seating

R82<sup>®</sup>



User height  
700–1800 mm



Max. user mass  
80 kg



Environment  
Indoor & outdoor

# Mustang<sup>®</sup>

Mustang is a highly adaptable walking aid for children and youngsters and it can be configured with anterior or posterior support to meet a wide range of walking needs.


The Mustang is available in 4 sizes.

## Highlights about Mustang



### Motivating the child to walk

The user is placed at the centre of the Mustang which provides core stability and enables the user to focus energy on walking.

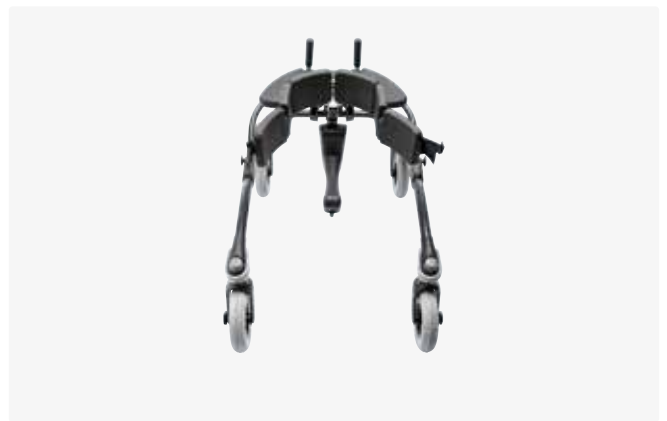

### Easy to access

The Mustang has been constructed with accessibility in mind. The U-shaped frame and wrap around support, allows for easy access and transfer.



### Natural walking pattern

The Mustang has an easy to adjust central bar for height and angle adjustment. A slight forward leaning position, combined with built-in suspension, stimulates the natural walking pattern.



### Accessories are easily added or removed for the carers convenience



## Standard

Four sizes, directional stabilizer, drag wheel function, tension adjustable wheels, angle adjustment, reverse brake on/off, height adjustment.

## Options

### Sizes

- Size 1
- Size 2
- Size 3
- Size 4

### Wheels

- 8" solid wheels
- 8" slime line wheels
- 4" wheels for size 1 + 2

### Frame colour



Yellow size 1



White size 2



Red size 3



Anthracite  
Grey size 4










Standard Mustang size 4



Standard Mustang size 1

## Measurements

|   |                             | 1          | 2          | 3           | 4           |
|---|-----------------------------|------------|------------|-------------|-------------|
|  | Seat height above floor*    | 260–440 mm | 330–500 mm | 450–680 mm  | 620–840 mm  |
|  | User height: floor to axil* | 400–650 mm | 600–850 mm | 700–1050 mm | 870–1300 mm |
|  | Max. user mass              | 30 kg      | 40 kg      | 60 kg       | 80 kg       |
|  | Length*                     | 820 mm     | 900 mm     | 1000 mm     | 1100 mm     |
|  | Overall width               | 520 mm     | 580 mm     | 650 mm      | 700 mm      |
|  | Minimum turning radius      | 460 mm     | 550 mm     | 600 mm      | 650 mm      |
|  | Mass device*                | 13 kg      | 14 kg      | 18 kg       | 19 kg       |

Further measurements at our website. \*Measured with 8" wheels. The Mustang is measured as standard delivered and with the center bar in vertical position. Different settings may result in another set of measurements.. To accomplish some of the measurements above, a slight adjustment might be required. We supply various types of accessories, which is why the weight is listed without

## Nice to know

Mustang is a highly adaptable walking aid, which can be used as a posterior or anterior walker. A recent systematic review was unable to make recommendations for the use of anterior or posterior walkers based on the existing evidence. However, the majority of walking aid user and their parents preferred posterior walkers (Poole 2017).

For children who are “new” walkers, the anterior support makes them feel safe and keep them from hitting objects with the body as they practice their walking skills. The more “experienced” walkers might feel free and more independent and able to get closer to objects with posterior support.




## Accessories



**Chest support**  
8698x0



**Cover**  
8698x0-x60

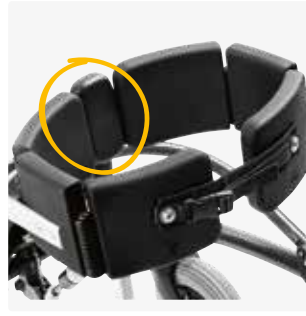

**Extension piece**  
869225



**Hand bar**  
86936x



**Abdominal support**  
869260



**Arm rests, anterior**  
86935x



**Arm prompts**  
8693xx




**Forearm positioner**  
869370



**Sacrum support**  
86925x



**Hip support**  
8692xx

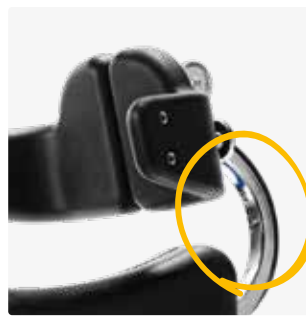

**Hip supports on curved bar**  
869222 or 869232



**Connecting rod**  
869226

For more accessories - visit our website



**Seat**  
8691x0



**Comfort seat**  
869115, 869125



**Sling seat**  
86957x



**Leg pad**  
869470



**Guide bar**  
869510



**Storage pouch**  
869560



**Storage bag**  
869565



**Leg separator**  
8694x0



**Ankle prompts**  
8694xx

